
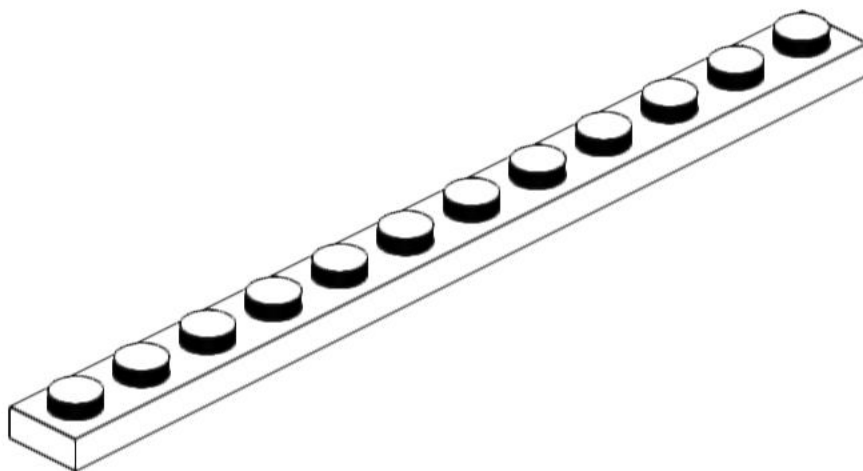




LEGOLAS Syringe Plunger Assembly Instructions

1

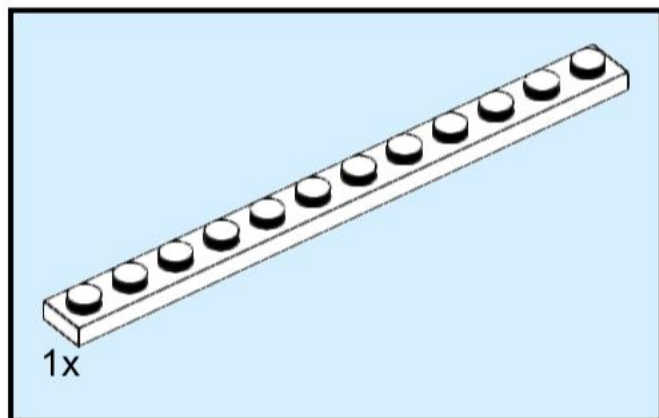



You may use a
1x12 plate LEGO
(as shown here)
or three 1x4
plates lined up

2

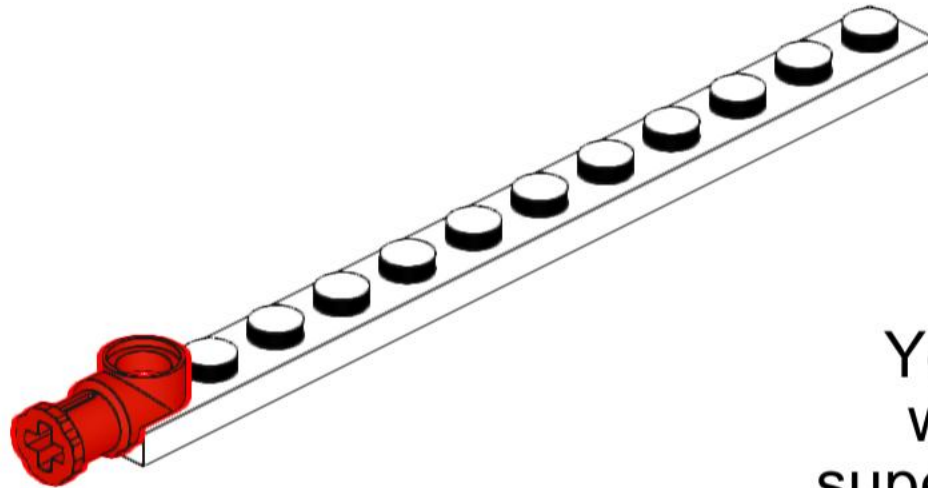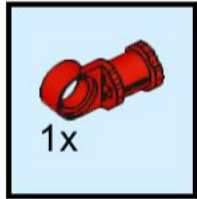


You may
want to
superglue the
toggle joint
onto the plate
button

3

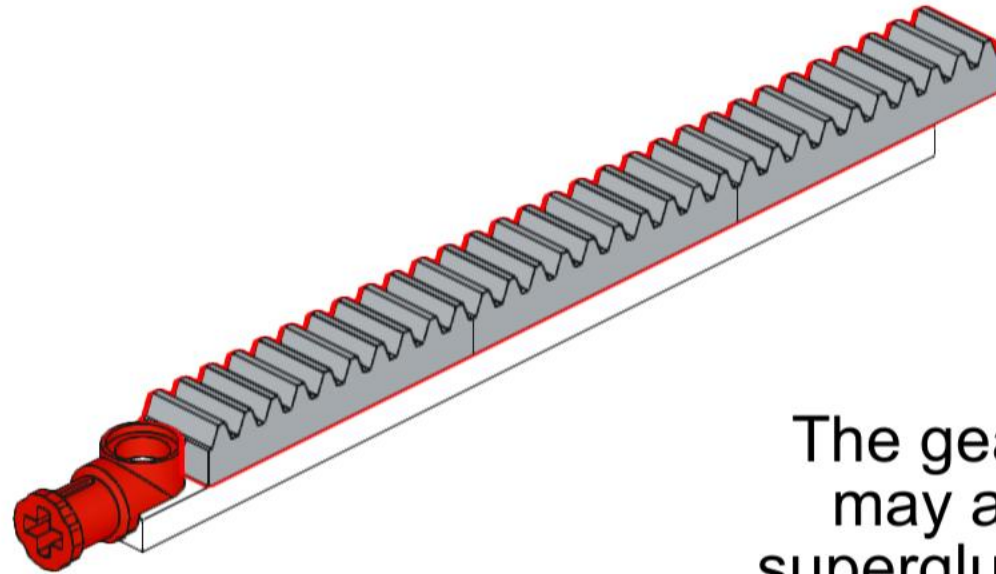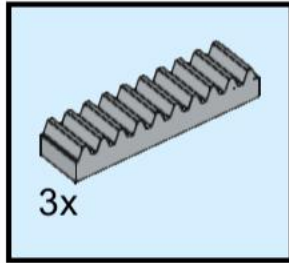

The gear racks
may also be
superglued to the
plate, although
this is not
necessary

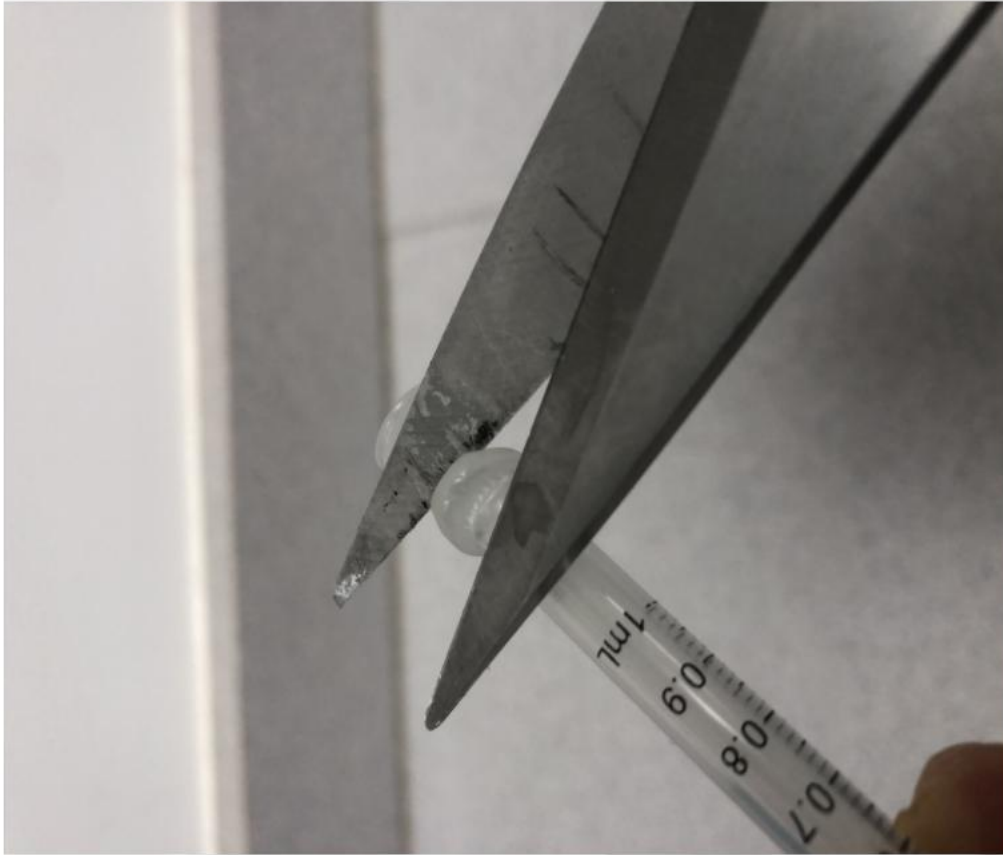

Use shears or scissors to cut off the syringe collar so that it has a thin profile (only a slightly larger diameter than the tube - as shown on the right)



Now cut off the plunger end (left), so that the 4-way cross is exposed, and insert this into one end of the toggle joint (right). Supergluing this joint is recommended.